





VITALITEC

cygnet®
Flexible Clamps

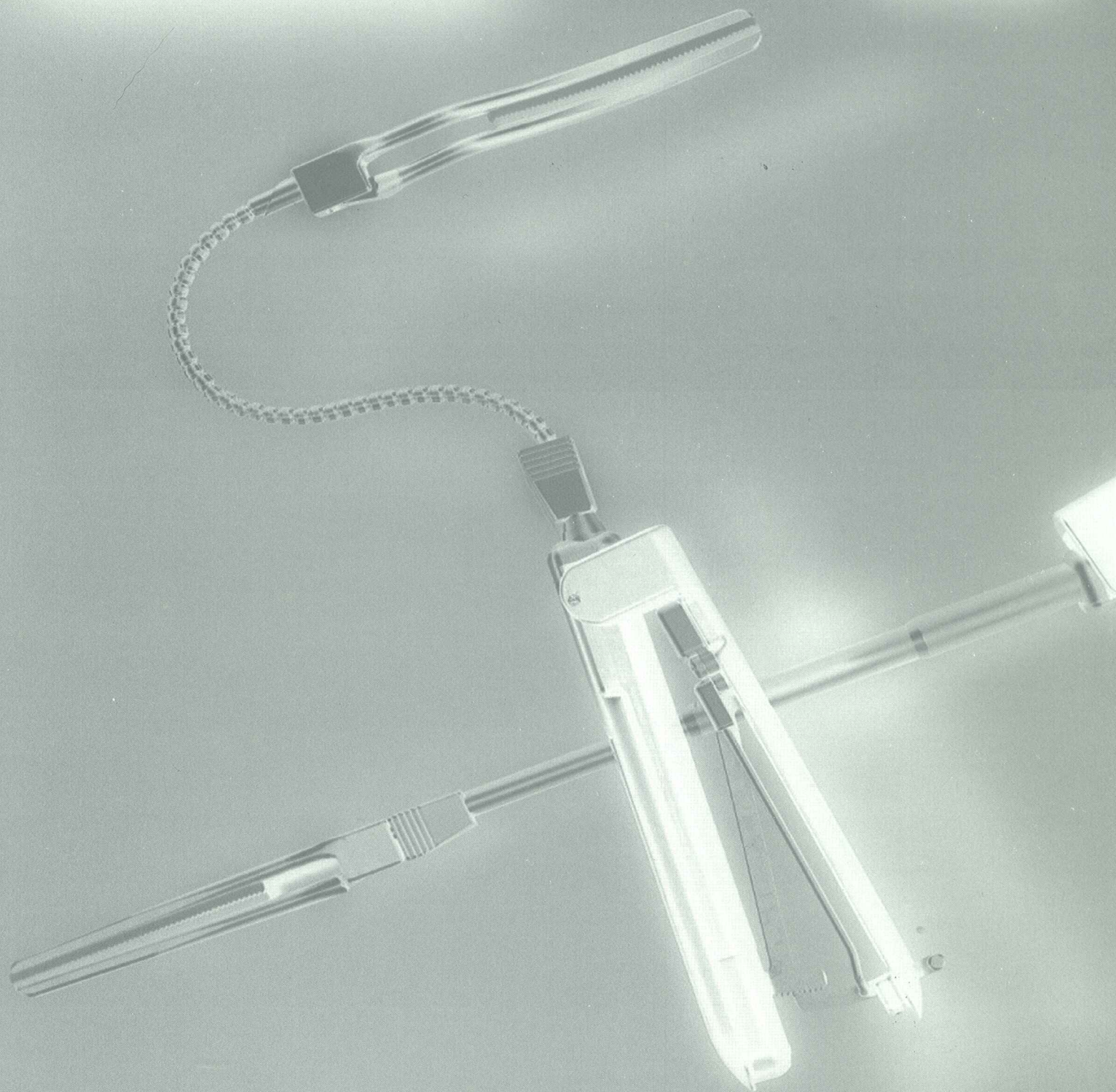

cygnet® captures the best of both worlds for multi-use and aortic cross clamping. The rigid shaft enables tunneling and accurate clamp placement. Once properly positioned, the rigid shaft retracts leaving a flexible neck, which can be maneuvered to a location remote of the surgical procedure.

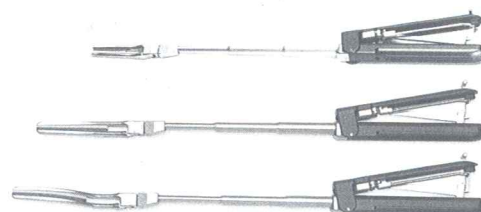


cygnet® Clamps

Name	Item #	Jaw Length	Jaw Orientation
Multi-Purpose Clamp	N-10183	33mm	Straight
Aortic Cross Clamp	N-10142	66mm	Straight
Arched Aorta Clamp	N-10143	66mm	Curve
Aortic Cross Clamp	N-10151	86mm	Straight
Arched Aorta Clamp	N-10152	86mm	Curve

intrack® inserts

	Item #	Length	Unit/Box
Conforming	N-10127	33mm	20/box
	N-10128	66mm	20/box
	N-10147	86mm	20/box
Standard/ Traction	N-10125	33mm	20/box
	N-10126	66mm	20/box
	N-10141	86mm	20/box
Double Traction	N-10121	33mm	20/box
	N-10122	66mm	20/box
	N-10146	86mm	20/box
Ultra	N-10174	66mm	20/box
	N-10175	86mm	20/box

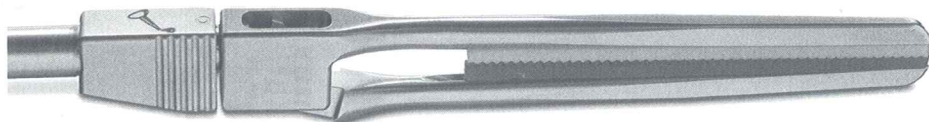

Controlled Strength

Comparable in strength to bulkier, aortic clamps, but able to provide greater access, control and flexibility in your procedures.

Inventive Design

The clamp jaw incorporates the patented intrack channel design, so you are assured your intrack inserts will remain secure in the channel and provide the highest level of vessel conformability and traction.

Ergonomic

An ergonomically designed handle gives you quick release capabilities and precision instrument ratcheting.

Orders within the U.S. may be placed by calling
Customer Service toll free 800.805.6590

re-order #:cyg0506-01

Vitalitec International, Inc.
10 Cordage Park Circle, Suite 100
Plymouth, MA 02360
508.747.6033 Tel
508.747.5118 Fax
www.VitalitecUSA.com
Email: VitalitecUSA@aol.com

Partner für Deutschland

FUMEDICA

Medizintechnik GmbH

Ziegelwasen 4
De-72336 Balingen
www.fumedica.de

Tel.: + 49 (0)7471 920 00 50
Fax.: + 49 (0)7471 920 00 52
fumedica@fumedica.de