



TITANIUM HEMOSTATIC CLIPS

CLIP 9 VITALITEC®



**Quality and Reliability
in one clip**

peters-surgical.com

Titanium hemostatic clips **CLIP 9 VITALITEC®**

The reference clip

Pure titanium


- Biologically inert.
- Annealed to a consistent malleability.
- No shape memory for occlusive closure.
- Radio opaque.
- Non-magnetic for use with MRI or related scanning equipment.



The cartridge

Ease of loading


- Adhesive backing holds the cartridge securely on any surface.
- Soft loading of the clip due to the unique design of the cartridge.
- Number of clips per cartridge better suited for number of clips used per procedure.

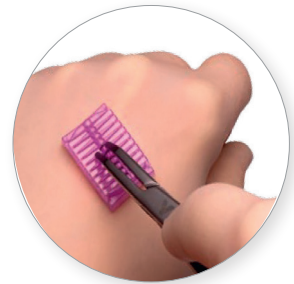

6 sizes, 1 exclusive size, 7 colour codes

- Clips sizes adapted to perform hemostasis on various vessel configurations.

9 Clips/Cartridge for Micro thru Medium and 6 Clips/Cartridge for Medium/Large and Large

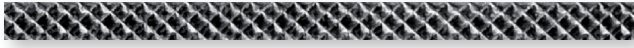


Internal structure

- Internal structure made up of small diamond shape with pyramidal heads.
- Interlocking atraumatic teeth designed for non-slip closure on the vessels.



Triangular shape wire

- The triangular shape of the titanium wire improves the precision of the positioning of the clip into the jaws of the applicator.

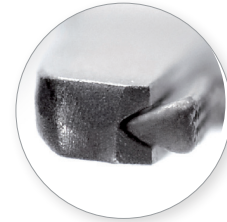


Appliers: the perfect fit

- Open surgery range: Lengths of 15 cm, 20 cm, 28 cm and 28 cm angled models.
- Laparoscopic range: Lengths 33cm and 45cm. Ø 5mm, 10mm and 12mm.



- Maximum opening of the applicators for a full cleaning.

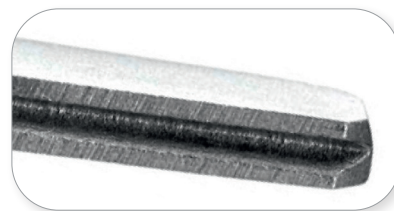

- Specific treatment of the inner surface of the jaws to secure the grip of the clip.











- Colour codes to fit all 7 available clips sizes.

Titanium hemostatic clips

CLIP 9 VITALITEC®



Titanium Hemostatic Clip 9 VITALITEC® | Cartridges






Size	Reference	Clips/Cartridge	Cartridge/Box	Clips/Box	Height clip closed 
 Micro	W9060	9	20	180	2.6 mm
 Small	J9180	9	20	180	3.6 mm
 Small	R9180	9	20	180	3.6 mm
 Small/Medium	L9180	9	20	180	4.7 mm
 Medium	B9180	9	20	180	5.6 mm
 Medium/Large	V9120	6	20	120	9.0 mm
 Large	O9120	6	20	120	12.3 mm

VITALITEC® | Hemostatic Clip Appliers

Size	15 cm	20 cm	28 cm	28 cm angled
<i>Jaws shape</i>				
 Micro	PW615	PW620	-	-
 Small	PJ115	PJ120	-	-
 Small	PR115	PR120	-	-
 Small/Medium	PL515	PL520	-	-
 Medium	PB215	PB220	PB228	PBC228
 Medium/Large	-	PV320	PV328	PVC328
 Large	-	PO420	PO428	POC428

References for open surgery appliers end with -EB, -E, -1B. References available depending on the country.
For more information please contact us.

VITALITEC® | Endoscopic Clip Appliers

Size	Reference	Required trocar Ø	Length
 Small	LCASLS1-1	Type 5 mm trocar	33 cm
 Small/Medium	LCASLS2-1	Type 5 mm trocar	33 cm
 Medium	LCASLS4-1	Type 10 mm trocar	33 cm
 Medium/Large	LCASLS5-1	Type 10 mm trocar	33 cm
 Large	LCASLS6-1	Type 12 mm trocar	33 cm

Other available references for endoscopic appliers, please contact us.

○..... The Titanium Hemostatic Clips 9 VITALITEC® are included in Peters Surgical products line○

Titanium Hemostatic CLIP 9 VITALITEC®

INDICATIONS: Titanium Hemostatic Clips can be used in every surgical specialty when complete occlusion of a vessel or tissue is required, except for Fallopian Tube Ligation.
Medical device class III, CE 0459. Manufacturer: PETERS SURGICAL

VITALITEC® Hemostatic Clip appliers

INDICATIONS: VITALITEC® Hemostatic Clip appliers are intended for the loading, holding and optimal closure of VITALITEC® Titanium Hemostatic clips.
Medical device Class I. Manufacturer: PETERS SURGICAL

VITALITEC® Endoscopic Clip appliers

INDICATIONS: VITALITEC® Endoscopic Clip appliers are intended for the loading, holding and optimal closure of VITALITEC® Titanium Hemostatic clips.
Medical device class I. Manufacturer: Vitalitec International Medizintechnik GmbH Germany distributed by PETERS SURGICAL

Presentation dedicated to PETERS SURGICAL employees, distributors and healthcare professionals.

Carefully read the Instructions For Use of medical device before use.



Peters Surgical Corporate headquarters: 42 rue Benoît Frachon 93013 Bobigny Cedex France | Tél. : +33 (0)1 48 10 62 62 | Fax : +33 (0)1 48 91 22 99
peters@peters-surgical.com - www.peters-surgical.com | RCS Bobigny 444 018 477



ENCLIP03v00D