
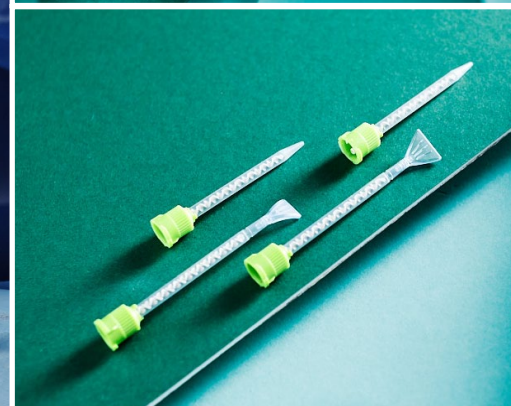


NEXGlue[®]

Surgical Adhesive

Your Safety Net

Full strength in **2 minutes**

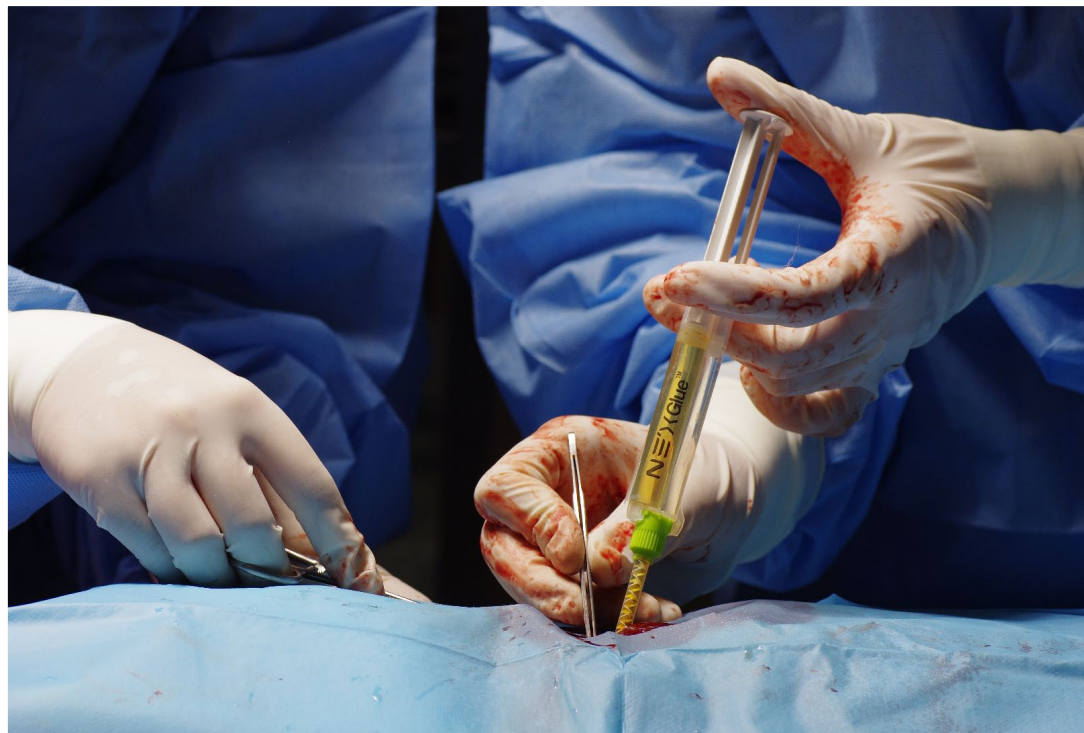

PROVIDING SPECIALISTS WITH
INNOVATIVE
HEALTHCARE SOLUTIONS

WITH THE CENTRAL FOCUS ON
PATIENTS' SAFETY AND USERS' COMFORT



NE'X Glue[®] Surgical Adhesive

NE'X Glue[®] is two component surgical adhesive which covalently bonds with tissue and takes a form of flexible hydrogel. It also mechanically binds with synthetic materials. Such features allow for strong tissue - tissue or tissue - synthetic material connections.

Combination of Bovine Serum Albumin (BSA) and glutaraldehyde in optimal ratio with precise mixing inside the applicator of the delivery device allows for full consumption of glutaraldehyde during the polymerisation. As a result, absolutely biocompatible flexible hydrogel is created.





KEY FEATURES:

- Variety of available volumes – 2, 5 and 10 ml
- Standard and spreader applicator tips
- Room temperature storage
- Ready for use in seconds
- Easy to use
- Full polymerization in 20-30 seconds
- Full strength in 2 minutes
- Works independently of the coagulopathic state of the patient
- Stronger than sealants based on fibrin or polyethylene glycols


CLINICAL APPLICATIONS:

NE'X Glue[®] is used as an adjunct to standard methods of repair such as sutures, staples, electrocautery and patches to seal, bond and reinforce tissue. Can be applied alone to fix surgical meshes, reinforce, repair or seal damaged parenchyma as a replacement of ineffective or impractical techniques.

- **VASCULAR SURGERY**
- **CARDIAC SURGERY**
- **THORACIC SURGERY**
- **GENERAL SURGERY**
- **GENITOURINARY SURGERY**
- **NEUROSURGERY**
- **ENT SURGERY**
- **SPINAL SURGERY**

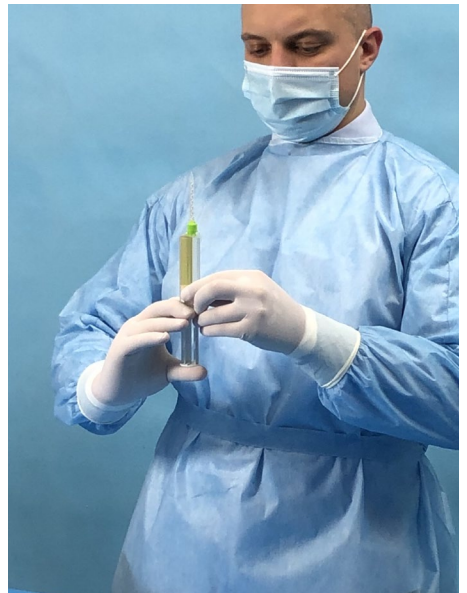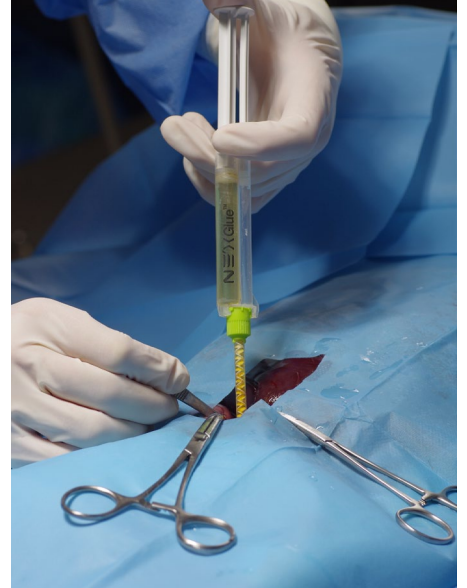

REF	PRODUCT	CONTAINS
0206-NX2	NE'X Glue® Surgical Adhesive 2 ml	5 single packs - each contains: 1 syringe (2 ml) & syringe plunger, 4 syringe applicator tips
0206-NX5	NE'X Glue® Surgical Adhesive 5 ml	5 single packs - each contains: 1 syringe (5 ml) & syringe plunger, 4 syringe applicator tips
0206-NX10	NE'X Glue® Surgical Adhesive 10 ml	5 single packs - each contains: 1 syringe (10 ml) & syringe plunger, 4 syringe applicator tips, 3 spreader tips (12 mm)
0206-NX4SM	Syringe applicator tip	10 single packs - each contains 4 syringe applicator tips
0206-NX3WM12	Syringe 12 mm spreader tip	10 single packs - each contains 3 syringe 12 mm spreader tips
0206-NX3WM16	Syringe 16 mm spreader tip	10 single packs - each contains 3 syringe 16 mm spreader tips

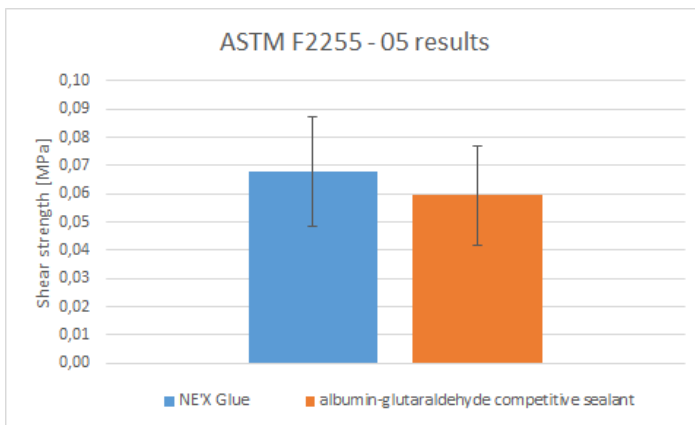

Endotoxin content in the whole extract [EU]	NE'X Glue® Surgical Adhesive	Albumin-glutaraldehyde competitive sealant
	0,290	0,535

Grena Biomed Ltd.
Chelsea House, Chelsea Street,
Nottingham, NG7 7HP, United
Kingdom

**CUSTOMER SERVICE CENTER
NOTTINGHAM OFFICE, UK**
☎ + 44 115 9704 800
☎ + 44 115 9704 800
✉ export@grena.co.uk
www.grena.co.uk

